Instruction Manual	2
Manuel d'instructions.....	8
Manual de Instrucciones	14
Bedienungsanleitung	21

Copyright 2005.
All rights reserved.
Printed in China.

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Safety Precautions

Sprout™ is a personal mini-centrifuge designed for laboratories in the bioscience, medical, and chemistry fields. It may be involved with biohazardous and/or hazardous materials. This instruction manual can not address all safety hazards. It is the responsibility of the user to consult and observe all health and safety precautions and to assess the instrument's suitability to the task.

Caution!

Before starting Sprout™ for the first time, please read this instruction manual carefully. This manual contains important information, safeguards, and operating instructions.

1. Never use Sprout™ in any manner inconsistent with these instructions.
2. Observe general laboratory safety precautions and regulations.
3. Do not attempt to operate with the cover removed or open.
4. The rotor and lid must always be securely fastened.
5. Always push the rotor all the way onto the shaft, then tighten the set screw.
6. Do not use damaged rotors.
7. Do not put hands in the rotor area unless the rotor is completely stopped.
8. Do not move unit while the rotor is spinning.
9. Do not immerse unit in water. Injury due to shock or fire may result.
10. The rotor must be loaded symmetrically. Operating this unit with an unbalanced rotor will cause severe vibration and damage, which is not covered by warranty.
11. Do not fill tubes while inserted in rotor. Spilled liquid may harm unit.
12. Close tube lids before starting centrifuge. Open lids can be ripped off during centrifugation and damage the unit.
13. Do not insert non-standard tubes or other foreign articles into the rotor.
14. Follow all precautions from the Material Safety Data Sheet (MSDS) for any reagent you use with this instrument.
15. Do not use solvent or flammable liquids near this or other electrical equipment.
16. Do not centrifuge flammable, explosive, or corrosive materials.
17. Do not operate in a hazardous or flammable environment.
18. Do not use any solvent on the unit that may attack plastic or cause cracks in the rotor.
19. Ensure the rotor is protected from corrosion and mechanical damage. The rotor must be cleaned with a pH-neutral cleaning liquid.
20. Always work in a manner which endangers neither the user nor any other person.

21. Should the instrument fail to work to specification, immediately stop centrifuging. Clean and troubleshoot the instrument according to the instructions under "Trouble Shooting" before any further use of the instrument. Contact an authorized dealer or the manufacturer if trouble persists.
22. Do not stop rotor while unit is turned ON. Doing so may cause the unit to fail and will void the warranty.
23. Only use original manufacturer's accessories (rotors, allen wrench, and adapters) and spare parts.
24. Repairs are to be performed by trained and authorized service personnel only. Contact Heathrow Scientific® LLC Repair Dept Technical Support.
25. Opening the instrument housing or improper use of the instrument voids the warranty. If there is a failure during the warranty period, contact Heathrow Scientific® LLC for warranty service.
26. Do not use units that have not been correctly installed or repaired.
27. Do not attempt to disassemble or modify.
28. Store at room temperature in a dry area. Do not expose to sunlight, moisture, or extreme temperatures for prolonged periods of time.

Technical Data

Maximum speed:	6,000 rpm
Maximum RCF:	2,000 x g
Electrical:	110/120 VAC @ 50/60 Hz or 220/230 VAC @ 50/60 Hz
Fuse:	0.10 A, 110 VAC

This mark is the confirmation that the unit conforms to the EU guidelines and has been tested according to the Low Voltage Directive 73/23/EEC.

Conforms to UL Std. 61010-1.
Certified to CAN/CSA C22.2 No. 1010.2.020-94.

Disposal Responsibilities:

The equipment you purchased may contain hazardous substances that could impact the environment. Per regulations on electronic devices in the European Community, you must use the appropriate disposal systems to avoid exposure of these substances to the environment. The disposal systems will reuse or recycle hazardous materials from your equipment responsibly.

The crossed-out wheeled bin symbol invites you to use those systems.

If you need more information on the collection, reuse and recycling systems, please contact your local or regional waste administration. You can also contact us for more information on the environmental performance of our products.

Operation

1. Remove the contents from the package, examining them carefully for breakage, defects, or missing parts. Save the packaging.
2. Place Sprout™ on a sturdy, level surface away from the edge of the counter and other moving equipment. Make sure the power cord is clear of hot surfaces and other hazards.
3. Verify that the power switch is in the OFF position. Plug the power cord into the rear of the unit and an approved, grounded outlet.
4. Turn power switch ON. If there is a smooth whirring sound and the unit accelerates with little or no vibration, Sprout™ is ready to use. If there are loud or unusual sounds, or excessive vibration, DO NOT OPERATE. Contact Heathrow Scientific® LLC Technical Support.
5. On benchtop prepare samples tubes in foam microtube rack provided so that tubes are filled to equal levels. Close tube lids.
6. To open lid, firmly tilt lid toward rear of unit.


Warning!

Spin balanced loads only!

Do not place the tubes asymmetrically nor load unequal volume tubes. Improper placement of tubes will lead to insufficient centrifugation and may cause injury or accident.

7. Tubes of equal weight and size should be placed opposite each other. Use sample or water in tube to provide balance. Proper sample balancing illustrated below.

One-Size Circular Rotor

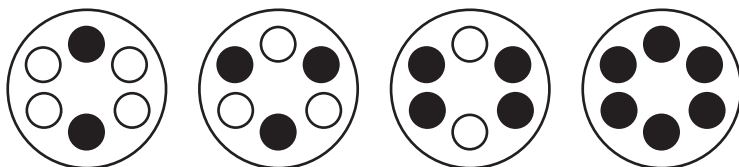
2 Tubes

3 Tubes

4 Tubes

6 Tubes

Dual-Size Circular Rotor

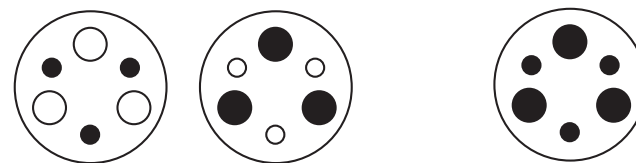
3 Tubes

3 Tubes

6 Tubes

Strip Tube Rotor

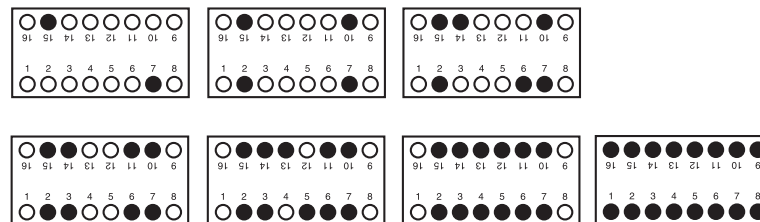8. Do not force tubes into the rotor. Use adapters and correctly sized rotors for your application.
9. Close lid securely. Do not attempt to operate with the lid removed or open.
10. Turn power switch ON. Unit will not operate with lid open. If there is a smooth whirring sound and the unit accelerates with little or no vibration, Sprout™ is ready to use. If there are loud or unusual sounds, or excessive vibration, DO NOT OPERATE. Contact the Heathrow Scientific® LLC Technical Support.
11. To stop centrifugation, press power switch and wait until the rotor has completely stopped rotating before opening the lid.

Warning!

Do not open lid until rotor has completely stopped!

Avoid severe personal injury or property damage from moving parts.

Rotor Replacement

1. Use included allen wrench to loosen set screw on rotor.
2. Pull rotor off.
3. Push new rotor all the way onto the shaft.
4. Tighten set screw.

Warning!

Do not use damaged or cracked rotors!

Maintenance

Sprout™ is normally maintenance-free. Clean unit only when it is not plugged into an electrical outlet. When necessary, the housing and rotor can be wiped using a damp cloth and a mild, non-corrosive detergent.

Warning!

Do not use any solvent on the unit that may attack plastic or cause cracks in the rotor. Ensure the rotor is protected from corrosion and mechanical damage. The rotor must be cleaned with a neutral cleaning liquid.

Notes

—Ensure all parts are thoroughly dry prior to operation.

Return for Repair

Important:

Transporting hazardous materials without a permit is a violation of federal law.

Heathrow Scientific® LLC will not accept any Sprout™ that is not appropriately cleaned and decontaminated.

In the unlikely event of repair, or when damage to the unit necessitates return, contact Heathrow Scientific® LLC and obtain return authorization **before** sending your Sprout™ for service.

Warranty

- ACCESSORY PARTS, SUCH AS RUBBER PARTS, THAT ARE DAMAGED BY LIQUIDS OR MISUSE.

IN NO EVENT WILL HEATHROW SCIENTIFIC® LLC'S OBLIGATION UNDER THIS WARRANTY EXCEED THE PRICE OF THE PRODUCT.

Limited Warranty

Heathrow Scientific® LLC warrants that Sprout™ will be free from defects in workmanship and material for 1 year from the date of purchase.

If you believe that there is a defect in the product, you must, during the warranty period, notify Heathrow Scientific® LLC, provide proof of purchase, and return the product to Heathrow Scientific® LLC with a Return Authorization form. To obtain a Return Authorization form, please call 1-847-816-5070. If Heathrow Scientific® LLC is properly notified and, after inspection, confirms that there is a defect and the warranty period has not expired, Heathrow Scientific® LLC will repair, modify, or replace the product, at its sole option, at no charge.

OTHER THAN THIS LIMITED WARRANTY, HEATHROW SCIENTIFIC® LLC MAKES NO WARRANTIES, EXPRESS OR IMPLIED, WITH RESPECT TO THE QUALITY OR PERFORMANCE OF THE PRODUCT, INCLUDING BUT NOT LIMITED TO THE IMPLIED WARRANTIES OF MERCHANTABILITY AND FITNESS FOR PARTICULAR PURPOSE WHICH ARE HEREBY DISCLAIMED AND EXCLUDED. HEATHROW SCIENTIFIC® LLC WILL IN NO EVENT BE LIABLE FOR ANY LOSS OF USE, LOSS OF PROFITS, CONSEQUENTIAL, SPECIAL, EXEMPLARY OR PUNITIVE DAMAGES.

THIS WARRANTY DOES NOT COVER:

- ANY DEFECT OR DAMAGE CAUSED BY IMPROPER OR UNREASONABLE USE OF THE PRODUCT. (SPROUT™ IS DESIGNED FOR USE ONLY BY TRAINED LABORATORY TECHNICIANS. USE BY ANYONE ELSE WILL VOID THIS WARRANTY.)
- ANY PRODUCT THAT HAS BEEN, IN HEATHROW SCIENTIFIC® LLC'S SOLE JUDGEMENT, TAMPERED WITH, ALTERED, OR REPAIRED BY ANYONE OTHER THAN HEATHROW SCIENTIFIC® LLC.
- ANY PRODUCT THAT IS INOPERATIVE BECAUSE OF: (a) WEAR OCCASIONED BY USE, (b) NEGLIGENCE, (c) ACCIDENT, (d) INCORRECT MAINTENANCE, OR (e) USE UNDER ABNORMAL CONDITIONS OF TEMPERATURE, DIRT OR CORROSION, OR USE WITH ABRASIVE OR CORROSIVE MATERIALS.

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